

SPECIFICATION

产品规格书

REFOND P/N 产品型号

RF-MTQAB7N03-B1

□R&D 研发

■Mass Product 量产供货

Contents 目录

1. Description 产品介绍	3
1.1 General Description 产品描述	3
1.2 Features 产品特征	3
1.3 Application 产品应用	3
1.4 Features Introduction 产品功能	4
2. Product Specification / 产品规格	5
2.1 Size chart 产品尺寸图	5
2.1 Product Parameters 产品参数	6
3 .application 成品应用	6
4 .Packing Criterion 包装标准	7
5. Handling Precautions 产品使用注意事项	8
5.1 Product storage and maintenance 产品存储和维护	8
5.2 Handling Precautions 产品使用注意事项	9

The logo for REFOND is centered on the page. It features a stylized graphic of several overlapping, curved, light blue lines above the word "REFOND" in a bold, sans-serif font. The entire logo is rendered in a light blue, semi-transparent watermark style.

REFOND

1. Description 产品介绍



1.1 General Description 产品描述

Notes 备注: The product consists of LED module (red light, blue light, IR), mask and control panel
产品由 LED 模组 (红光、蓝光、IR) 和面罩以及控制板组成的 LED 面部美容仪

1.2 Features 产品特征

- ▶ Appearance of private module design, fashion and beautiful;
外形私模组设计，时尚美观;
- ▶ According to the facial uniform design, so that every inch of skin feel light;
根据面部均匀设计，使每一寸肌肤感受到光;
- ▶ Low voltage DC3.7V constant voltage driving to ensure safe and reliable running ;
低压 3.7V 工作，确保运行安全可靠；
- ▶ Visual design, beauty while working;
可视设计，边工作边美容;
- ▶ Built-in battery, easy to carry ;
内置电池，方便携带.
- ▶ Touch switch, switch at will.
触控开关，随心切换.

1.3 Application 产品应用

- ▶ Facial beauty app

1.4 Features Introduction 产品功能

▶ Red light, 650-660nm, promotes blood circulation, accelerates metabolism, promotes collagen production under the skin, and reduces fine lines and wrinkles;

红光 650-660nm，促进血液循环，加速新陈代谢，促进皮下胶原蛋白的生成，减少细纹和皱纹；

▶ Blue light 405-420nm, improve the pore environment, improve acne, reduce the frequency of acne;

蓝光 405-420nm，改善毛孔内环境，改善痘痘，减少发痘频率；


▶ Powder light +IR (red light & Blue light & IR) can improve the alternating function of cell oxygen, improve microenvironment, fade pigment, soften wrinkles and moisturize skin ;

粉光+IR (红光&蓝光&IR) 改善细胞氧的交替功能，改善微环境，淡化色素，抚淡皱纹弹润肌肤.



2. Product Specification / 产品规格

2.1 Size chart 产品尺寸图

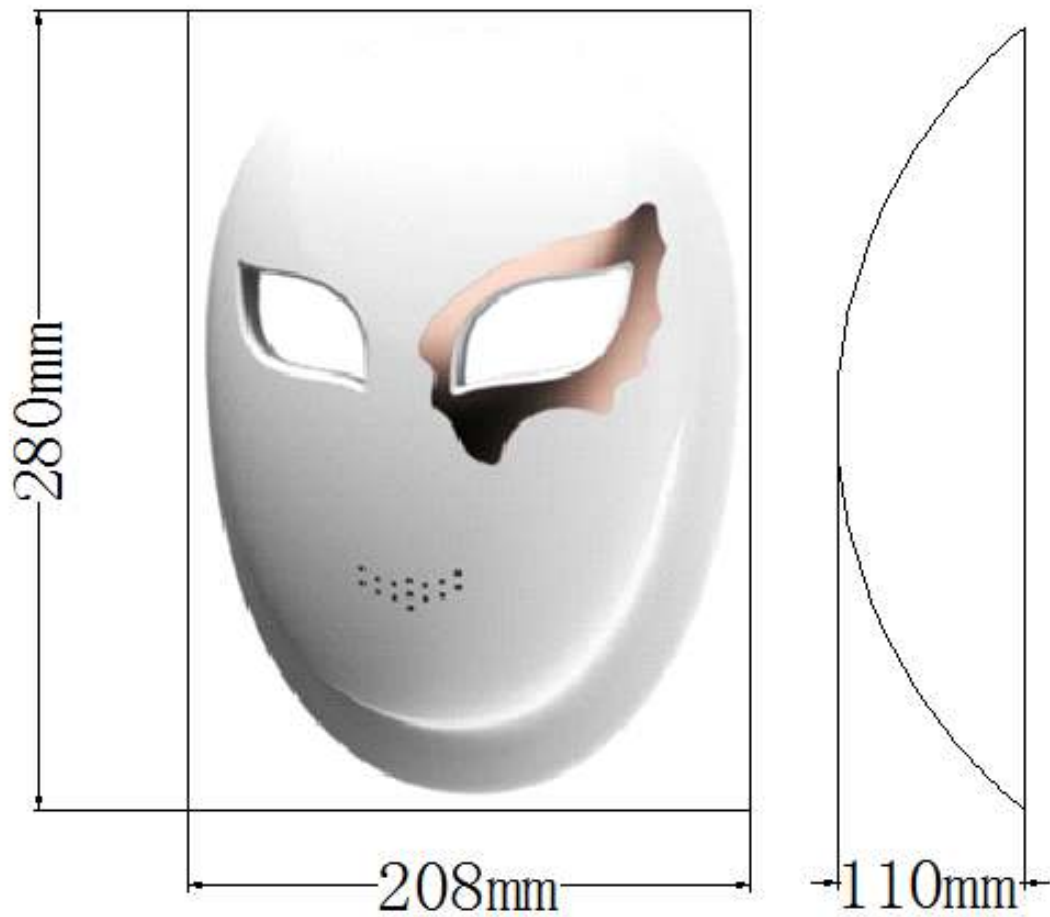

尺寸size : 280*208*110mm

2.1 Product Parameters 产品参数

产品参数 (Product Parameters)		
序号	名称	产品参数
1	Refond PN 瑞丰产品料号	RF-MTQAB7N03-B1
2	Battery 电池	3.7V/2300mAh
3	Touch 触控	模式一：粉光 (RED+IR+UV) 点亮 20 分钟 pink light (red + IR + UV) is on for 20 minutes
		模式二：蓝光 (UV) 点亮 20 分钟 blue light (UV) is on for 20 minutes
		模式三：红光 (RED) 点亮 20 分钟 red light (red) is on for 20 minutes
4	LED	RED:650-665nm
		UV: 405-420nm
		IR: 850-860nm
5	LED 数量 Number of LEDs	203PCS
6	充电器 Charger	开关电源适配器 电压 5V 电流 2A Charger switching power adapter Voltage 5V current 2A
7	充电头 Charging head	TYP-C
8	额定功率 Power	8W
9	净重 Net weight	
10	产品包装尺寸 Product package size	330*250*140mm
11	整箱尺寸 Full carton size	680*520*300mm
12	整箱数量 full carton quantity	8PCS

3 application 成品应用

模式一：粉光 Mode 1: pink light



模式二：蓝光 Mode 2: Blu light



模式三：红光 Mode 3: red light




4 Packing Criterion 包装标准

4.1 Package Diagram 包装示意图



Fig.3-2 Title 图名

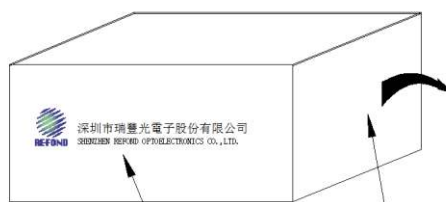

正反面印上“瑞丰LOGO”等图案

两侧面印上“料号、规格”和
纸箱三防标志等图案

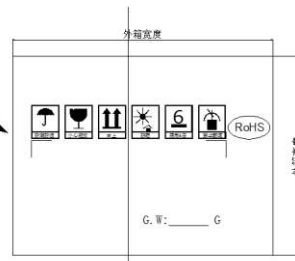


Fig.3-2 Title 图名

5. Handling Precautions 产品使用注意事项

5.1 Product storage and maintenance 产品存储和维护

▶ Do not expose the product to rain, humidity and high temperature.

请勿将本产品暴露在雨中、潮湿、高温环境中。

▶ Open flame sources, such as lit candles, should not be placed on the device

装置上不应放置明火火源，如点燃的蜡烛。

▶ Careful supervision is required to prevent children from using this product.

需要小心监督以防止儿童使用本产品。

▶ If not used for a long time, charge the battery at least once every 3 weeks to maintain the battery life

如果长时间不使用，请至少每 3 周对电池充电一次以维持电池寿命。

▶ Do not poke this product with sharp metal objects such as keys, hairpins, necklaces or nails.

请勿用尖锐的金属物体（如钥匙、发夹、项链或指甲）戳本产品。

▶ Do not put this product in the oven, dishwasher, washing machine

请勿将本产品放入烤箱、洗碗机、洗衣机中。

5.2 Handling Precautions 产品使用注意事项

Products for facial care, not recommended for body use.

▶ 产品用于面部护理，不建议用于身体。

Do not use this device for skin allergies, wounds, severe acne, burns, and bacterial skin infections.

▶ 皮肤过敏、创伤或严重粉刺、灼伤、细菌性皮肤感染等，请勿使用本设备

Please consult a dermatologist before using this device if you are being treated for a skin disease or have a skin condition.

▶ 正在接受皮肤病相关治疗或有皮肤问题，请先咨询皮肤科医生，然后再使用本设备。

Not recommended for pregnant, menstruating, lactating women (unstable hormonal balance may cause temporary skin problems).

▶ 不推荐怀孕、月经、哺乳期妇女使用(不稳定的荷尔蒙平衡可能会引起暂时性的皮肤问题)。

Do not use for people who are unable to express themselves clearly (including babies).

▶ 请勿无法清楚表达自己意志（包括婴儿）的人使用。

People with light sensitivity should not use this product if they are allergic to ultraviolet light.

▶ 有光敏感的人如紫外线过敏，则请勿使用该产品

Do not use this product in the rain, damp (bathtub, shower, sink), high temperature (open flame source, oven) environment.

▶ 请勿在雨中、潮湿（浴缸、淋浴器、水池）、高温（明火源、烤箱）环境中使用本产品。

Do not disassemble without permission, there is no user repairable parts inside, related maintenance, please consult qualified maintenance personnel.

▶ 请勿私自拆卸，内部无使用者可维修部件，相关维修，请咨询合格的维修人员。

Version History/修订历史

Date日期	Revisor修订者	Version版本	Verifier审核	Remarks备注
2021-8-11	范锐文	E/0	邱涛	The first edition

